


ImplantBook®

The ultimate global guide 2018



DENTAL IMPLANTS

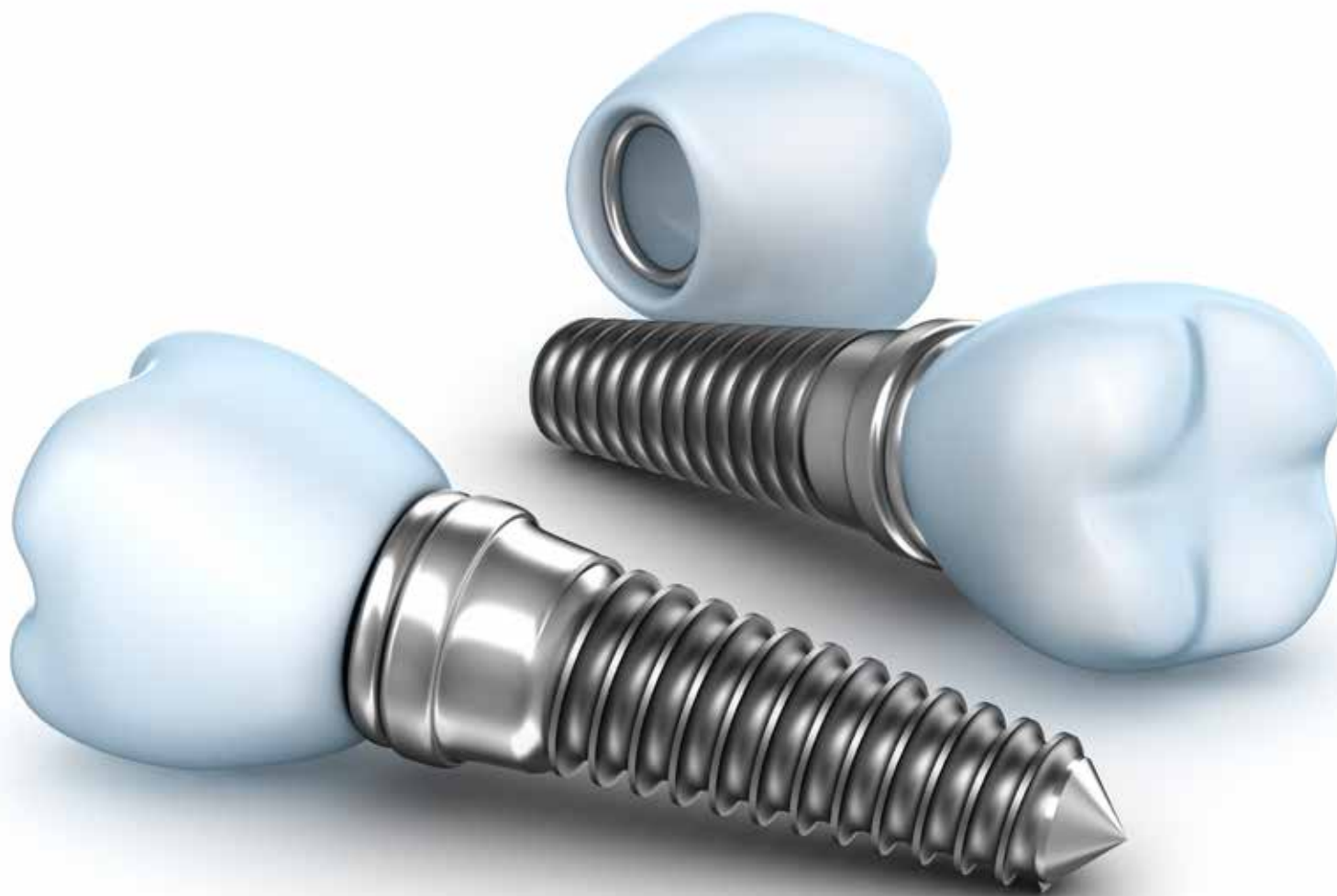

OSSEOINTEGRATION
& BIOMATERIALS



AROUND
IMPLANTOLOGY



3D PRINTERS

EQUIPMENT

PIEZOSURGERY

SOFTWARE

RADIOLOGY

CONGRESSES



COMPANY LIST



• Dental Tech

ImpLogic® GII



Manufacturer: ImpLogic® GII, evolution of the previous ImpLogic®, produced by Dental Tech, completes the safe and wide range of products that the company makes available to meet the needs of all professionals. The unique prosthetic platform unites component engagement, facilitating prosthetic procedures and reducing management costs.

Implant Morphology: ImpLogic® GII is a conical implant system particularly indicated in cases where it is necessary to safeguard structures anatomic as the maxillary sinus and alveolar nerve. The geometry of the spire aids osseous healing, both qualitatively and quantitatively. The micro-thread of the collar reduces stress while screwing the implant in the cortical area of the bone and contributes to the increase in primary stability and to the maintenance of bone level.

Surface: The BWS® surface, made by Dental Tech through a sandblasting and etching process, provides an extremely clean surface and an optimum roughness that promotes cell adhesion.

Available Implant Lengths (mm): 8.5 (only for \varnothing 4.5), 10, 11.5, 13 and 14.5 (only for \varnothing 3.75). The polished crown portion range between 0.3 (\varnothing 4.5 and \varnothing 5.5) and 1 (\varnothing 3.75).

Available Implant Diameters (mm): \varnothing 3.75, \varnothing 4.5 and \varnothing 5.5.

Seating Position To The Bone Crest: It's a submerged implant with three helical discharge grooves form three progressive cutting areas, providing the implant a marked ability to penetrate, while simultaneously removing bone debris, thereby reducing bone compression.

Abutment-Implant Connection Type: Hexagonal internal connection.

Index Type Abutment Repositioning: Hexagonal

www.dentaltechitalia.com // info@dental-tech.it